

S.no. Volume of NaOH added (ml)

Conductance ( $\Omega^{-1}$ )

01	0	5.6
02	0.5	5.0
03	1.0	4.3
04	1.5	3.5
05	2.0	2.9
06	2.5	2.4
07	3.0	2.2
08	3.5	2.4
09	4.0	2.5
10	4.5	2.7
11	5.0	2.8
12	5.5	3.0
13	6.0	3.2
14	6.5	3.6
15	7.0	4.0
16	7.5	4.5
17	8.0	4.9
18	8.5	5.4
19	9.0	5.8
20	9.5	6.3
21	10.0	6.7
22	10.5	7.1
23	11.0	7.5
24	11.5	8.0
25	12.0	8.4
26	12.5	8.8
27	13.0	9.2
28	13.5	9.6
29	14.0	10.0
30	14.5	10.4
31	15.0	10.7




Expt 7


Aim: To estimate the strength of mixture of acetic acid and hydrochloric acid by conductivity

### Principle:

Conductivity of a ~~soln~~ solution is related to size of and mobility of ion present when a mixture of strong acid (HCl) and weak acid ( $\text{CH}_3\text{COOH}$ ) is titrated against a strong base (NaOH), HCl reacts first followed by  $\text{CH}_3\text{COOH}$ , when titration of strong acid and strong base is carried out, there is a decrease in conductivity as highly mobile  $\text{H}^+$  ions are replaced by  $\text{Na}^+$ .

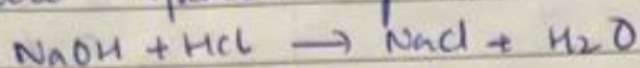

When the whole strong acid is consumed, base reacts with weak acid and conductivity increases as unionized weak acid becomes ionized salt.

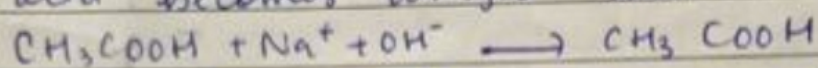

After both the acids are consumed, there is a sharp increase in conductivity which gives end point and this increase in conductivity is due to fast moving  $\text{OH}^-$  ion from NaOH solution.

By the amount of base consumed, amount of acid present is calculated.

### Procedure:

The solution is diluted to 100ml. 20 ml of which is pipetted out into a clean beaker and 80ml of distilled water is added. Conductivity cell is dipped into the test solution and titrated against NaOH with stirring. Conductance is measured for each



addition of 0.5 ml NaOH.

After neutralization, amount of acid present is determined by amount of NaOH consumption for strong acid and weak acid, volume of base consumed for strong acid and weak acid are determined by plotting a graph between conductance and volume of base added. First end point corresponds to strong acid while the other is for weak acid.

Calculation.

Strength of HCl

Volume of mix = 20 ml

Normality of HCl =  $P = N_1$

Vol. of NaOH = 2.7 ml =  $V_1$

normality of NaOH = 0.5 N

$$\text{Strength of HCl} = \frac{V_1 \times 0.5}{20} = 0.0675 \text{ N.}$$

Strength of  $\text{CH}_3\text{COOH}$

volume of mixture = 20 ml

normality of  $\text{CH}_3\text{COOH}$  =  $? = N_1$

Volume of NaOH = 3.3 ml =  $V_2 - V_1$

normality of NaOH = 0.5 N

$$\text{Strength of } \text{CH}_3\text{COOH} = \frac{0.5 \times (V_2 - V_1)}{20} = \frac{0.5 \times 3.3}{20}$$

$$= 0.0825 \text{ N.}$$

Result:

The strength of HCl present in given soln = 0.0675 N

The strength of  $\text{CH}_3\text{COOH}$  present in given soln = 0.0825 N



Exp 7

Scale  
X axis = 1 unit  
Y axis = 1 unit

